

Installation Instructions

Telescopic Arm for Concept LR

PROCEDURE:

1. Remove the existing set screw from the cuspidor arm.
2. Disconnect the cuspidor umbilical from the "J" box.



NOTE: Be sure to note the location of the umbilical tubing connections inside the "J" box. After the installation, the umbilical tubing will have to be reconnected in the "J" box.

3. Remove the cuspidor and umbilical from the cuspidor arm by lifting the cuspidor up and pulling the umbilical out of the arm (See Figure 1).

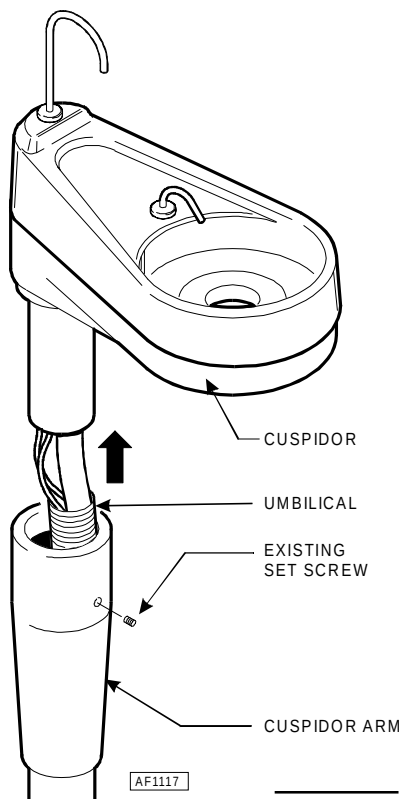


Figure 1

4. Position the support ring on the cuspidor arm and install the new adjusting set screw. Make sure the set screw goes in the hole where the old set screw was located (See Figure 2).
5. Level the support ring and tighten the two locking set screws.
6. Assemble the remaining parts and squeeze them together.
7. Tighten the three set screws in the cover ring.
8. Test the movement of the assembly and loosen or tighten the new adjusting set screw if necessary.
9. Replace the cuspidor and umbilical and reconnect the umbilical in the "J" box.

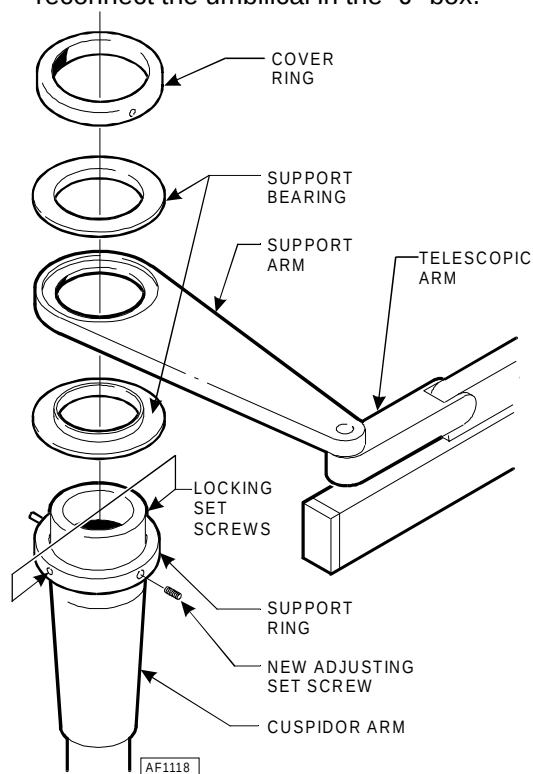


Figure 2

Knight®
by Midmark

Midmark Corporation • P O Box 286 • 60 Vista Drive
Versailles • Ohio • 45380-0286 • www.midmark.com
Phone: 937-526-3662 Fax: 937-526-5542

TPW1723 Rev. B